


BILLY BOY CLASSIC Advertising Condom BILLY BOY®

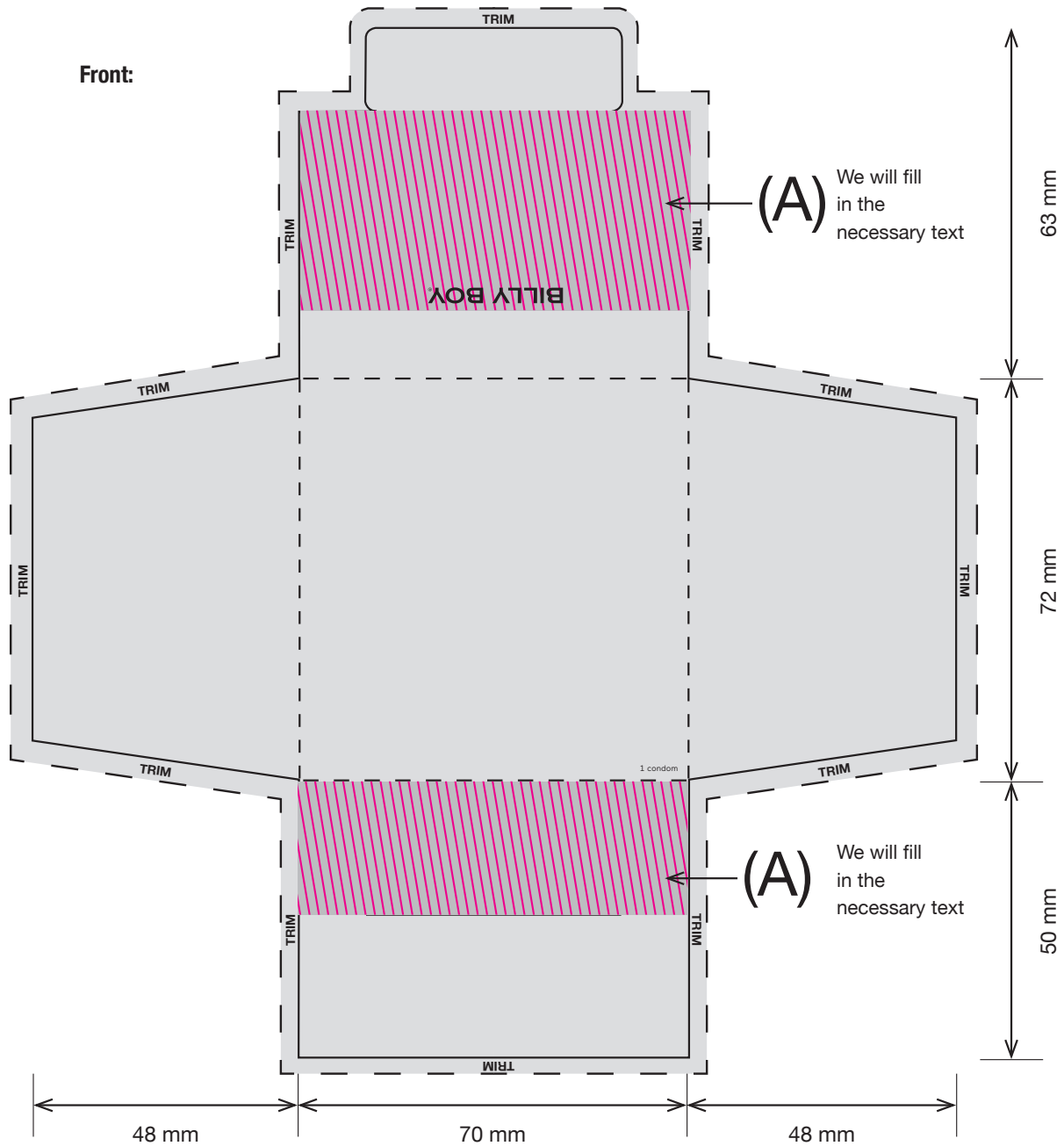

GUIDANCE

(A) = hatched areas in red:

Please don't fill in any data in these areas, they may not be written-on, changed or moved. We will insert the correct mandatory information. With exception of the BILLY BOY - Logo, the text colour can be chosen as per request.

Empty spaces:

All free spaces as well as the background colour can be designed individually based on CMYK Euro-scale (please consider the bleed as per marks). Special colours are only available on request. The element „1 condom Made in Germany for Sex“ can be moved to another corner of the front. Please remove all contour- and folding-lines!

Printing data:

PDF X-3 or -4, JPG, TIFF or EPS – format, vector files, min. 300 dpi resolution. The colour application in the file should not exceed 230 %.

Important notes:

MAPA GmbH does not assume any responsibility for legal infringements resulting from the use of material, protected by copyright or otherwise, which has been provided by the customer for the design of the advertising condom packs or which are used at his explicit request. We therefore recommend to check all design-elements correspondingly before usage. We reserve the right to change the layout or text of the mandatory information at any time, if necessary with regard to the identification of medical devices. For printing-approval, you will receive a non-colour-binding PDF - file with all mandatory condom data, implemented from our side. The printing data will be checked against our company policy guidelines and principles and we reserve the right to reject printing layouts even without giving reasons.











BILLY BOY CLASSIC Advertising Condom BILLY BOY®

Back side (inner side) - must not be used for your design!

All mandatory information and instructions for use for medical devices, which are required on the packaging by law, will be printed by us in black (on the white background) in the respective valid version and respective national language. You will find the explanation to this topic in our guidelines for placing condoms into the market. If you intend to export the condoms, all mandatory information has to be issued in the respective, official national language. Please ask for the export conditions and availability of the necessary language in advance.

If changes of the inner layout of the packaging are requested - e.g. adding any background colour or additional elements to the layout, based on CMYK - this will cause additional printing charges and you have to ask for a separate offer.

Please note the following information symbols on the packaging:

-  Please read the instructions for use carefully.
-  Protect condoms from direct sunlight and heat.
-  Store condoms in a dry place.
-  Do not use the condom if the wrapper is damaged.
-  Each condom should only be used once.
-  Condoms made from natural rubber latex.
-  Medical Device (EU) 2017/745
-  Manufacturer
-  Lot number
-  Expiry date

Instructions for use

1. Opening
Carefully tear open the condom wrapper at the marking ("open here"). The wrapper can thus be torn open cleanly without damaging the condom. Caution: The condom may be damaged by sharp objects, jewellery or fingernails. Do not use the condom if the wrapper is damaged or already open.

2. Before intercourse
Hold the condom test between the thumb and forefinger, making sure that no air is trapped inside. Pull back the foreskin and unroll the condom all the way down the erect penis before any contact of the penis with the body of your partner. If the condom does not unroll easily, it may be on backwards. In this case, do not turn it around for use as your partner may come into contact with body secretions. Throw it away and use a new condom, making sure that it is on the right way for unrolling. If you notice that the condom slips or tightens, check it for potential damage as this may cause a tear.

3. After intercourse
Withdraw the penis before it becomes flaccid! Make sure to hold the condom against the base of the penis to prevent the condom from slipping off.

Important notes:
Condoms should always be handled with care. Please check the expiry date on the package and the condom-wrappers - do not use the condom if expired! Store condoms in a cool and dry place and protect it from direct sunlight.
Condoms do not interfere with the human organism and have no side effects. There is, however, a possibility of intolerances occurring, which could involve allergic reactions (e.g. latex allergy) or even anaphylaxis. Should you encounter any discomfort or skin irritation, stop using the condom and consult your doctor. If used correctly (see instructions for use), BILLY BOY condoms provide protection from pregnancy and transmission of diseases. However, no contraceptives completely eliminate the risks of pregnancy or infection. If the condom has torn or broken during use, seek medical advice as soon as possible, but within 72 hours.

Please report any serious incident in relation to the device to the manufacturer and your national competent authority. Used condoms should be disposed of in a waste bin, not down the toilet. Wash penis and hands before resuming intimate contact. Condoms are designed to be used only once.

Attention:
Do not use condoms with oil-based or greasy lubricant (e.g. vaseline, baby-oil, bodylotion, massage-oils, cooking-oils). Those affect the condom and cause it to tear quickly. If using a lubricant, always ensure that the label states that it is suitable for use with condoms. Prior to using lubricants, please read the instructions for use on the packaging. Never apply lubricant directly onto the penis or the inside of the condom as this will increase the risk of the condom slipping off or tearing. When using condoms for anal intercourse, you should always also apply lubricant (e.g. BILLY BOY lubricant gel) to the condom as this reduces the risk of the condom slipping off or tearing. If vaginal suppositories are used, please comply with the respective instructions for use. If prescription or non-prescription medication (e.g. ointments, gels etc.) is applied to either the penis or vagina, please check with your doctor or pharmacist for compatibility of the medication with the condom. Caution: Condoms can be damaged by jewellery or body piercings. Each condom should only be used once. Reusing condoms makes them less safe and can lead to unwanted pregnancy and/or infection with HIV or contagious sexually transmitted diseases. Use a new condom in case of any change of sexual activity (vaginal, anal, oral).

